


The Grace ALTO System includes all the essentials needed to save time in the OR and greatly simplify middle ear surgery.



Disposable Sizers



ALTO Cutter

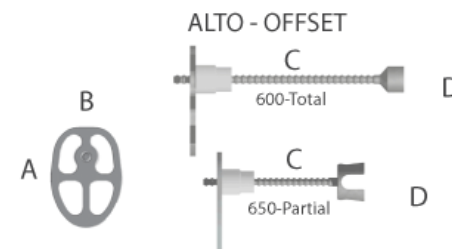

Cat. #	Material	A	B	C	D
600	Titanium	2.50mm	3.50mm	3.00-7.00mm	0.8mm
650	Titanium	2.50mm	3.50mm	2.00-5.60mm	1.45mm

The Grace ALTO line offers a unique patented method that makes the ALTO the only non-crimp adjustable length titanium design in otology today. This advantage allows the surgeon to adjust the prosthesis even after initial sizing, trimming and placement. The prosthesis can be lengthened or shortened in precise .25mm increments.