

OR/Mobius®

Maximize exposure, safety and efficiency with the OR/Mobius Self-Retaining Elastic Abdominal Retractor for Minilaparotomy




Improved minilaparotomy safety, efficiency and speed are just some of the benefits of the OR/Mobius® retractor

Once inserted, using a simple roll down technique, the Mobius® retractor provides 360° of radial retraction: maximizing working area and providing significantly greater exposure than traditional retractors.


The Mobius retractor comes in three sizes optimal for general, colorectal and gynecologic surgeries. Each size distributes retraction force evenly without excess pressure around the incision, creating unobstructed exposure to the surgical field. Throughout the procedure this device: maximizes working area; minimizes trauma, reduces risk of wound infection and blood loss, while promoting patient recovery.

OB/Mobius retractors are also available for use in c-section procedures. For more information please visit www.coopersurgical.com/obmobius

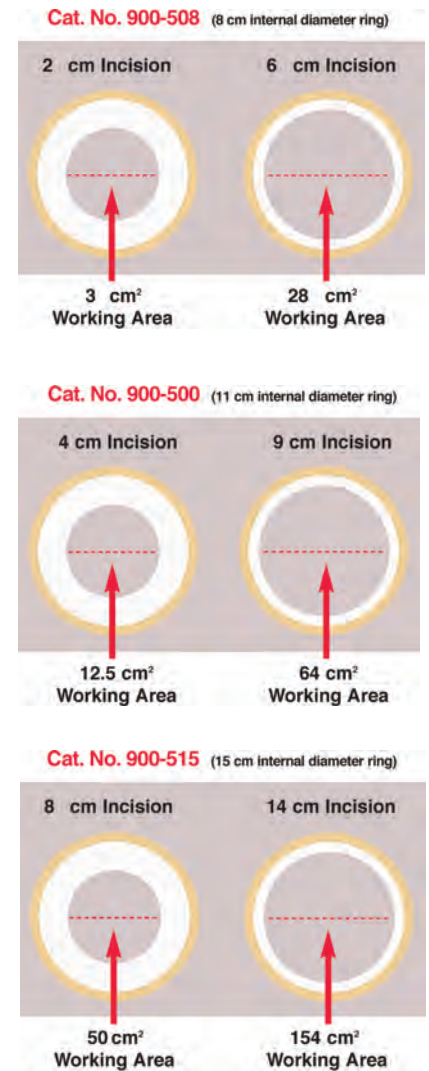

Simple Roll-down Technique:

1. Squeeze the lower blue ring and insert into the abdomen. Gently pull the ring to ensure insertion and perform a quick digital exam to check for trapped tissue.
2. Position hands at 11 and 1 o'clock with assistant at the 12 o'clock position.
3. Roll down the ring outwards first with the left, followed by the right, slide fingers down the ring to complete a Mobius flip. Repeat process 2-3 times until desired exposure achieved.

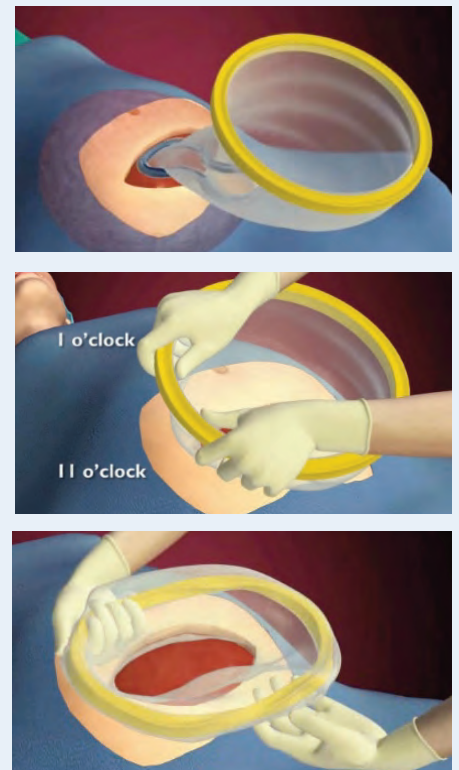

A better retractor means better visualization.


Our patented outer ring provides better flexibility and form.

Better Flexibility

- Outer ring lies flush against patient's skin without excessive pressure reducing pain and potentially reducing blood loss.
- Outer ring lies flush against patient's skin without interfering with surgeon's wrists.
- Flexible polyurethane sheath maintains moisture at the incision site while reducing risk of wound infection.

Better Form

- One-piece, molded outer ring construction rolls outward providing maximum hands-free exposure even in obese patients.
- Molds gently to the patient's abdomen for secure low profile retraction that won't slip.
- Only 2-3 outward roll-downs required to achieve optimal exposure.
- Tear resistant, non-stretch sleeve provides superior tissue retraction while allowing for movement through the wound.
- Bottom ring maintains shape, preserving exposure throughout procedure.

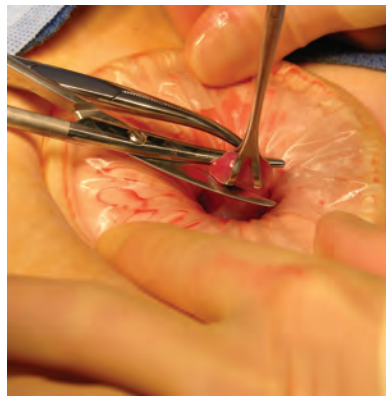

Easy to place.

Easy to use.

Easy to remove.

Benefits of the Mobius:

- Maximizes exposure through radial retraction for enhanced visibility
- Presents unobstructed surgical field unencumbered by traditional stainless steel retractors
- Alleviates need for additional surgical staff
- Protects wound edges minimizing contamination^{1,2}
- Gentle retraction through dispersed tension promotes patient postoperative comfort
- Ideal for obese patients (including BMI > 30)
- Easy to place and remove with simple technique

Applications:

- Cystectomy
- Total Abdominal Hysterectomy
- Myomectomy
- Tubal Ligation
- Colectomy
- Gastric Bypass
- Thoracotomy - Paraesophageal hiatal hernia
- Appendectomy
- Splenectomy

OR/Mobius

| PRODUCT NUMBER / DESCRIPTION: | | |
|-------------------------------|--|----------------------------------|
| 900-508 | Mobius, Mini inner diameter 8 cm | Optimal incision size 2-6 cm |
| 900-500 | Mobius, Small inner diameter 11 cm | Optimal incision size 4-9 cm |
| 900-515 | Mobius, Medium inner diameter 15 cm | Optimal incision size 8-14 cm |

Also available in larger sizes for maximizing exposure during c-sections.

OB/Mobius®

| PRODUCT NUMBER / DESCRIPTION: | | |
|-------------------------------|--|-----------------------------------|
| 900-517 | Mobius, Medium-Large inner diameter 17 cm | Optimal incision size 9-15 cm |
| 900-520 | Mobius, Large inner diameter 20 cm | Optimal incision size 12-17 cm |

Mobius is latex-free, made of 100% medical grade polyurethane. It is supplied sterile, one retractor per pouch. 5/box

References:

- ¹ Minimally Invasive Cesarean: Improving an Innovative Technique. Pelosi & Pelosi, OBG Management, July 2004
- ² The Use of the Mobius® Self-Retaining Retractor in Cesarean Deliveries—A Retrospective Case. Finck, MD & C Della Badia, Poster, AAGL 33rd Annual Meeting, Nov. 2004

Proudly distributed by



1300 00 ENDO (3636)

endotherapeutics.com.au

CooperSurgical

CooperSurgical is the leading company dedicated to providing medical devices and procedure solutions that improve health care delivery to women regardless of clinical setting. Our company is fostering that position through expansion of its core businesses and introduction of advanced technology-based products which aid clinicians in the management and treatment of commonly seen conditions.

203.601.5200 | 800.243.2974 | coopersurgical.com | 95 Corporate Dr. | Trumbull, CT 06611
also at www.YouTube.com/CooperSurgical

Copyright © 2011 CooperSurgical

Form #82019 Rev 07/11