

Danis Procedure Pack

SX-ELLA Stent Danis



Danis



Removable stent system for treatment of bleeding esophageal varices

Easy implantation without the need for endoscopy / X-Ray control

Secure hemostasis through direct compression of varices

Oral intake maintained from time of implantation



3D animation



Basic description

Procedure Pack with SX-ELLA Stent Danis

Self-expandable nitinol stent preloaded in a ready-to-use delivery system. Complete procedure pack including all items required for implantation.

Indications

SX-ELLA Stent Danis is intended to stop acute and/or refractory bleeding from esophageal varices.

Contraindications

- ☑ Diagnosed malignancy or stricture affecting the esophagus
- ☑ Case history of irradiation therapy of the chest or the esophagus
- ☑ Malignancy of the throat or larynx
- ☑ Malignancy of bronchus or fistula
- ☑ Stomach malignancy
- ☑ Suspicion of bleeding due to foreign body injuring the upper gastrointestinal tract

Features / Benefits

- ☑ Variable stent braiding is adjusted to esophageal peristalsis - reduced migration rates.
- ☑ Standardized compression of the varices results in effective hemostasis.
- ☑ Atraumatic stent ends.
- ☑ Radiopaque markers at both stent ends and at the midpoint - for possible X-ray control of the stent position.
- ☑ Retrieval loops made of medical grade alloy at both stent ends - high mechanical durability and acid resistance.
- ☑ Easy endoscopic stent removal within 7 days.

Benefits of the method

- ☑ Readily implantable without endoscopy / X-Ray control, also in emergency situations.
- ☑ Patient can continue an oral diet after stent implantation.

Additional description


Danis Procedure Pack (019-08S-25-135) consists of a plastic carry case containing all items required for stent implantation:

- ☑ Delivery system preloaded with stent, guide wire, plastic syringes, bite block, paper emesis basin, slobber cloth, waste bag with plastic clip.

Danis Procedure Pack Basic (019-08S-25-135-B) consists of:




- ☑ Delivery system preloaded with stent, guide wire, plastic syringe 50 ml.

The delivery system 9.4 mm (28 F) / 6.6 mm (20 F) is available in a standard length of 60 cm.

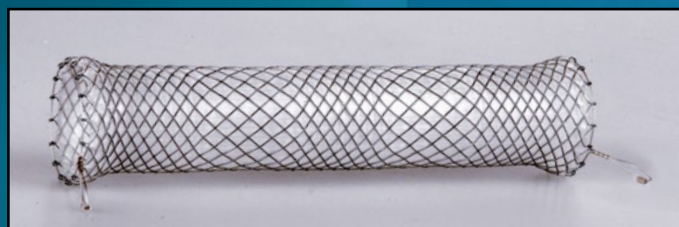

Available sizes

REF No.	SX-ELLA Stent Danis			Delivery system	
	Stent flares diameter [mm]	Stent body diameter [mm]	Nominal length [mm]	Active length [cm]	Outer diameter [F]
019-08S-25-135	30	25	135	60	28 / 20
019-08S-25-135-B	30	25	135	60	28 / 20

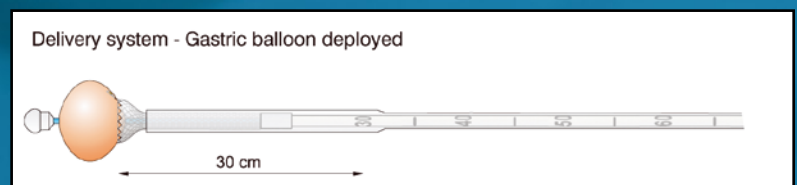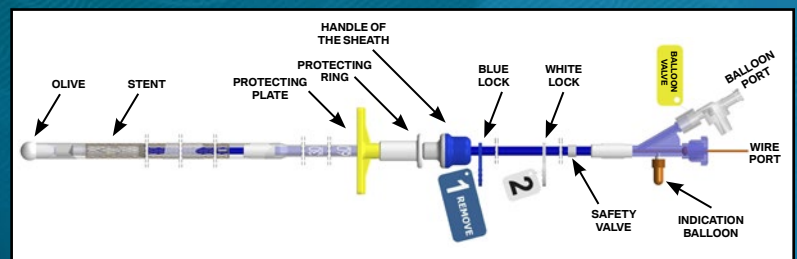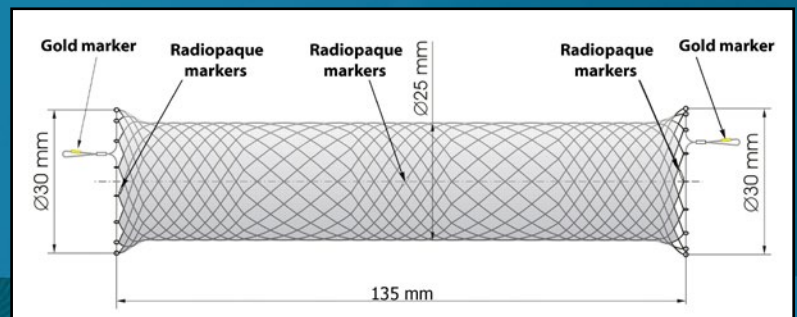

Proudly distributed by

Endotherapeutics
IMPROVING HEALTHCARE

☎ 1300 00 ENDO

✉ info@endotx.com.au

🌐 endotx.com.au

Manufacturer:

ELLA-CS, s.r.o.

Milady Horákové 504/45, Třebeš
500 06 Hradec Králové, Czech Republic

Tel.: +420 495 279 111 E-mail: info@ellacs.eu

URL: www.ellacs.eu