


# HV Stent Plus

SX-ELLA Stent Esophageal HV

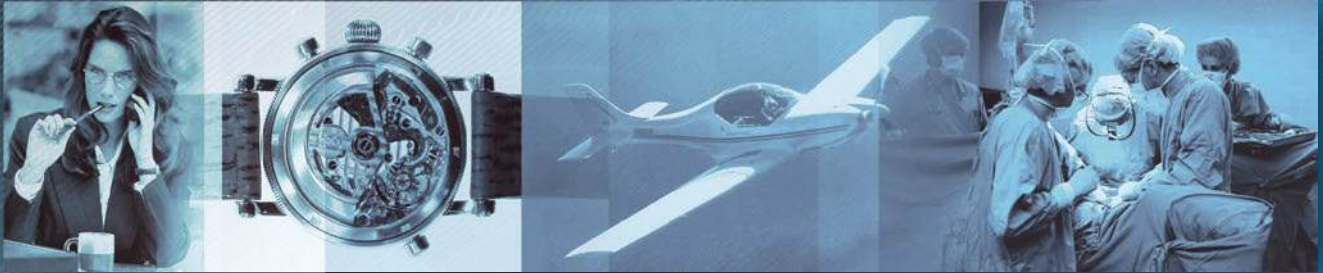

Durable silicone covering

Removable esophageal stent

Special anti-migration design

Atraumatic flares

Excellent flexibility



## Basic description

### SX-ELLA Stent Esophageal HV – HV Stent Plus

The HV Stent Plus is intended for patients aged 18 and older and indicated for palliative treatment of malignant dysphagia, for sealing malignant esophagorespiratory or esophagomediastinal fistula. The device can be temporarily implanted in patients with refractory benign esophageal strictures.

Stent system for fluoroscopic and endoscopic insertion.  
Self-expandable fully covered nitinol stent preloaded in PULL delivery system.




### Features / Benefits

- ☑ Retrieval loop made of medical grade alloy available on both stent ends - high mechanical durability and acid-resistance; stent flare centripetal constriction; suitable for extraction by ELLA Extractor.
- ☑ Radiopaque markers - accurate stent positioning.
- ☑ Anti-migration collar\* - reduced migration rates.
- ☑ Atraumatic stent ends - esophageal wall irritation and hyperplasia reduction.
- ☑ Excellent flexibility - smooth alignment and adequate stent expansion along anatomical curves.
- ☑ Durable silicone stent covering - tumor in-growth prevention, esophagorespiratory fistulas occlusion.
- ☑ Optional anti-reflux valve - gastroesophageal reflux and aspiration prevention.

\* Designed by: UMC Utrecht, The Netherlands

### Additional description

5.2 mm / 7.2 mm (16 F / 22 F) delivery system available in standard length of 75 cm.  
Recommended guide wire: 0.035" (0.89 mm) - ultra stiff - minimum length of 220 cm.

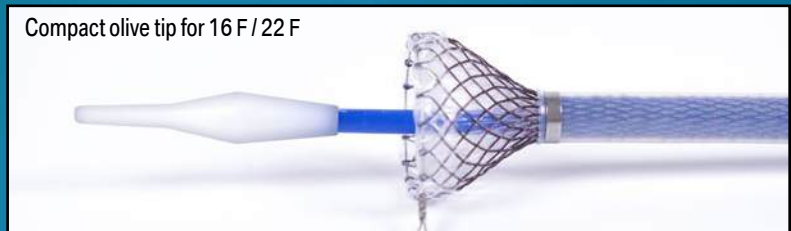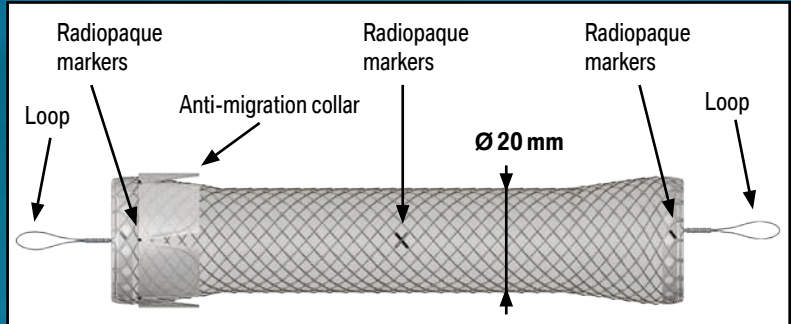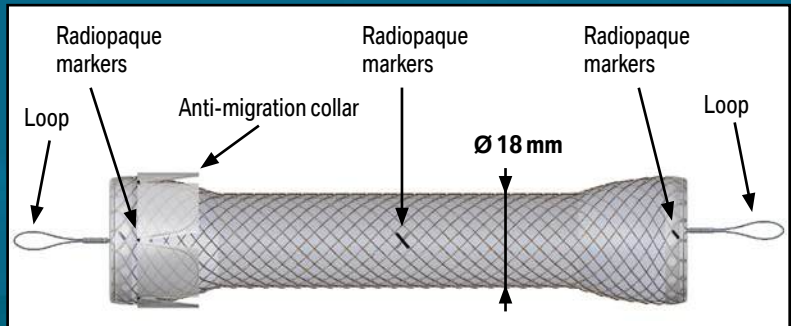

Delivery system - 16 F / 22 F



### Available sizes

REF No.	Stent diameter [mm]		Stent nominal length [mm]	Anti-reflux valve	Delivery system
	Proximal flare / Body / Distal flare				
019-09SL-20-085	25/20/25		85	No	16 F / 22 F
019-09SL-20-110	25/20/25		110	No	16 F / 22 F
019-09SL-20-135	25/20/25		135	No	16 F / 22 F
019-09SL-20-150	25/20/25		150	No	16 F / 22 F
019-09SL-20-085-O	25/20/25		85	Yes	16 F / 22 F
019-09SL-20-110-O	25/20/25		110	Yes	16 F / 22 F
019-09SL-20-135-O	25/20/25		135	Yes	16 F / 22 F
019-09SL-20-150-O	25/20/25		150	Yes	16 F / 22 F
019-09S-18-085	25/18/25		85	No	16 F / 22 F
019-09S-18-110	25/18/25		110	No	16 F / 22 F
019-09S-18-135	25/18/25		135	No	16 F / 22 F
019-09S-18-150	25/18/25		150	No	16 F / 22 F
019-09S-18-085-O	25/18/25		85	Yes	16 F / 22 F
019-09S-18-110-O	25/18/25		110	Yes	16 F / 22 F
019-09S-18-135-O	25/18/25		135	Yes	16 F / 22 F
019-09S-18-150-O	25/18/25		150	Yes	16 F / 22 F

Anti-migration collar



Proudly distributed by

**Endotherapeutics**  
IMPROVING HEALTHCARE

☎ 1300 00 ENDO  
✉ info@endotx.com.au  
🌐 endotx.com.au

Manufacturer:

**ELLA-CS, s.r.o.**

Milady Horákové 504/45, Třebeš  
500 06 Hradec Králové, Czech Republic

Tel.: +420 495 279 111 E-mail: info@ellacs.eu  
www.ellacs.eu