

Coloproctology
Urology
Urogynecology
Laparoscopy
Infusion systems
Obesity

i-Cut

Ergonomic & efficient morcellation

i-Cut is the user-friendly single-use morcellator with reliably high cutting performance. Ideal for extracting large tissue masses in LASH and myomectomy.


Convenient. Powerful. Economical.



Urogynecology

A.M.I.[®]

Agency for Medical Innovations

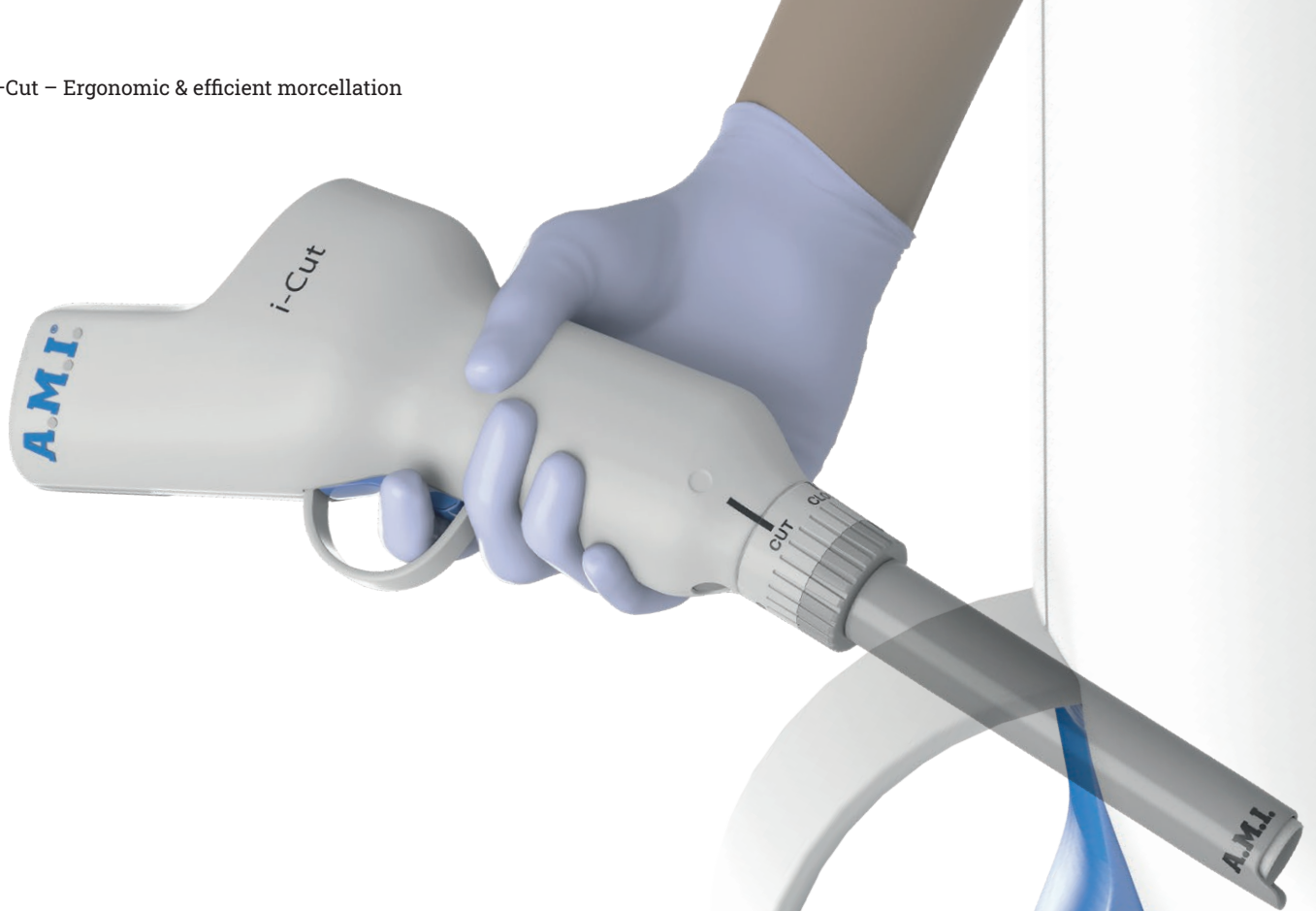

i-Cut

Developed by users for users.

The i-Cut disposable morcellator combines the best of both worlds: All the advantages of a maintenance-free single-use solution plus consistently high cutting performance even during operations on large and calcified tissue.

Technical data

Speed: 1000 rpm
Diameter: 15 mm
Shaft length: 100 mm
Weight: approx. 400 g
Power: 24 V DC, 160 W

User-friendly

- Ergonomically shaped handpiece
- Lightweight
- Plug-and-play: easy setup and handling in the OR

Powerful

- Precision ground blade
- No battery, therefore constant high cutting performance
- Large cutting tube diameter for the removal of large and hard tissue masses



More-Cell-Safe

More-Cell-Safe is the first solution available on the market that significantly reduces the risk of unexpected cancer cell scatter during laparoscopic morcellation. Optic and morcellator are introduced via two separate accesses. Small pieces of tissue, cell particles and fluids resulting from morcellation are isolated in a closed, cell-safe system. In addition, the special optics sleeve Visi-Shield reduces contamination of the optics. In combination with the powerful i-Cut, safety during surgery is further increased.



Working lumen

Activation button

Handpiece

Adjusting wheel

Trocar

Cutting tube

Made in Europe

Precision-ground cutting-tube.
Developed and manufactured in Austria.

Safe

- Multi-stage safety concept to protect the blade – “no-cut position” and safety clip to prevent unintentional exposure of the cutting tube
- Double trocar sealing system prevents gas loss

Economical

- Time and cost efficient one-step solution
- Fast operating
- No preparation costs for sterilization, maintenance and service
- Low operating costs

Order number	Product	Description
ICT5011 	i-Cut Power Morcellator for extracting tissue in LASH and myomektomy	i-Cut consists of: Power morcellator and obturator, 1 piece, delivered sterile
ICT5020 	i-Cut Power Supply Type C – Euro Reusable power supply unit for i-Cut ICT5011	1 piece, delivered non-sterile Additionally available: ICT5030 i-Cut Power Supply Type G ICT5040 i-Cut Power Supply Type I ICT5050 i-Cut Power Supply Type A ICT5060 i-Cut Power Supply Type C - Korea
MCS5111 	More-Cell-Safe 0° System for reducing the scattering of cancer cells during morcellation during laparoscopic surgery For 0° optics	The More-Cell-Safe system consists of: <ul style="list-style-type: none"> • Cell-tight, intraperitoneally insufflatable More-Cell Bag with two accesses • Visi-Shield to prevent camera contamination with 0° optics 10 systems per box, delivered sterile
MCS5151 	More-Cell-Safe 30° System for reducing the scattering of cancer cells during morcellation during laparoscopic surgery For 30° optics	The More-Cell-Safe system consists of: <ul style="list-style-type: none"> • Cell-tight, intraperitoneally insufflatable More-Cell Bag with two accesses • Visi-Shield to prevent camera contamination with 30° optics 10 systems per box, delivered sterile

The products contained in this brochure meet the requirements of the Medical Device Directive 93/42/EEC and are accordingly marked with the CE mark.

EN 200052

ICT5011, MSC5111, MSC5151

CE 0297

Proudly distributed by



☎ 1300 00 ENDO (3636)
 info@endotherapeutics.com.au
 www.endotherapeutics.com.au
 Suite 301, 2 Banfield Road
 Macquarie Park NSW 2113

A.M.I. Headquarters:

A.M.I. Agency for Medical
 Innovations GmbH
 Im Letten 1
 6800 Feldkirch
 Österreich
 T +43 5522 90505 0
 F +43 5522 90505 4006
 e info@ami.at

A.M.I.®

www.ami.at