

Coloproctology  
Urology  
**Urogynecology**  
Laparoscopy  
Infusion systems  
Obesity

# i-Cut

## Ergonomic & efficient morcellation

i-Cut is the user-friendly single-use morcellator with reliably high cutting performance. Ideal for extracting large tissue masses in LASH and myomectomy.


Convenient. Powerful. Economical.



Urogynecology

**A.M.I.**<sup>®</sup>

Agency for Medical Innovations

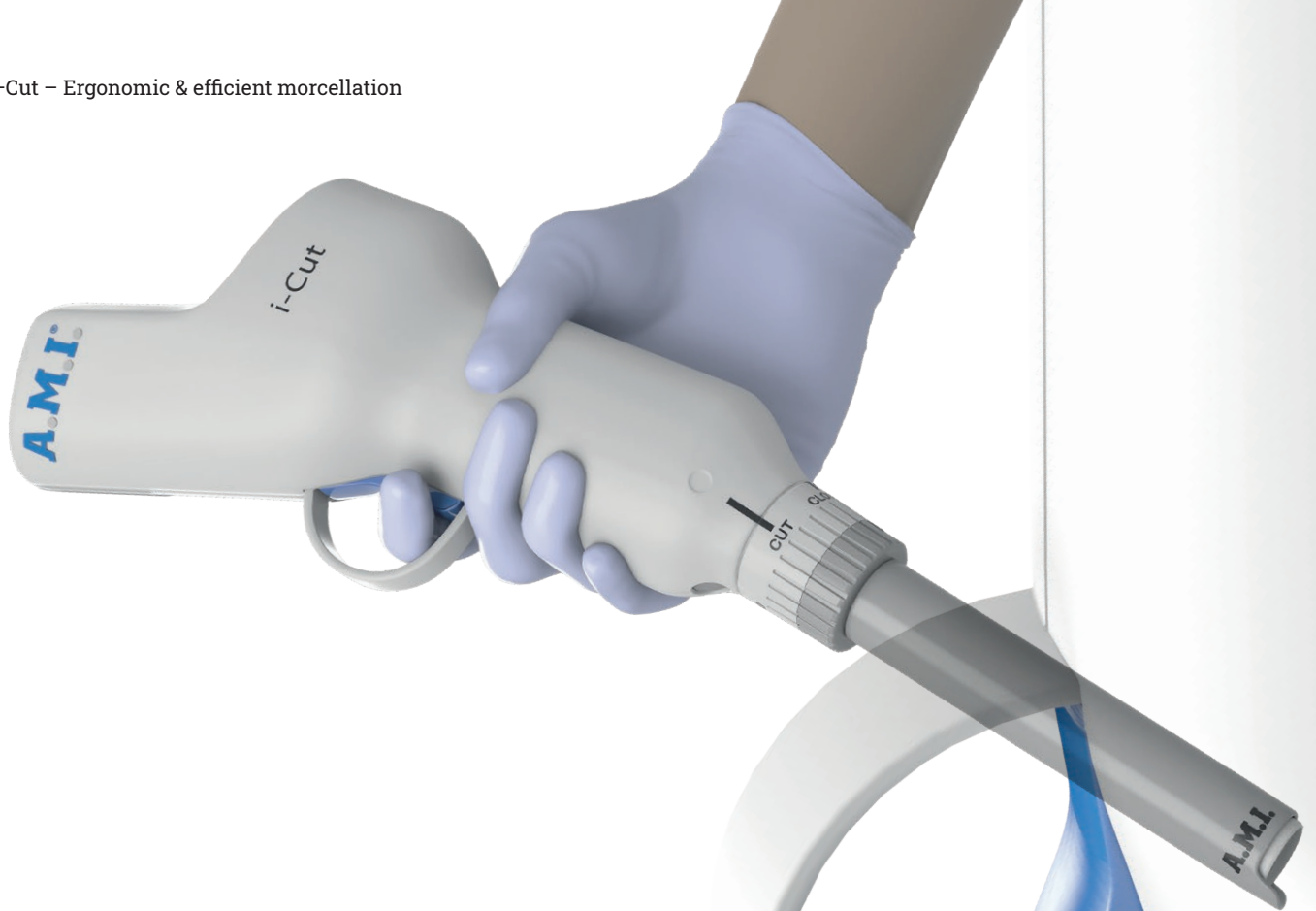

# i-Cut

Developed by users for users.

The i-Cut disposable morcellator combines the best of both worlds: All the advantages of a maintenance-free single-use solution plus consistently high cutting performance even during operations on large and calcified tissue.

## Technical data

Speed: 1000 rpm  
Diameter: 15 mm  
Shaft length: 100 mm  
Weight: approx. 400 g  
Power: 24 V DC, 160 W

## User-friendly

- Ergonomically shaped handpiece
- Lightweight
- Plug-and-play: easy setup and handling in the OR

## Powerful

- Precision ground blade
- No battery, therefore constant high cutting performance
- Large cutting tube diameter for the removal of large and hard tissue masses



### More-Cell-Safe

More-Cell-Safe is the first solution available on the market that significantly reduces the risk of unexpected cancer cell scatter during laparoscopic morcellation. Optic and morcellator are introduced via two separate accesses. Small pieces of tissue, cell particles and fluids resulting from morcellation are isolated in a closed, cell-safe system. In addition, the special optics sleeve Visi-Shield reduces contamination of the optics. In combination with the powerful i-Cut, safety during surgery is further increased.



Working lumen

Activation button

Handpiece

Adjusting wheel

Trocar

Cutting tube

## Made in Europe

Precision-ground cutting-tube.  
Developed and manufactured in Austria.

### Safe

- Multi-stage safety concept to protect the blade – “no-cut position” and safety clip to prevent unintentional exposure of the cutting tube
- Double trocar sealing system prevents gas loss

### Economical

- Time and cost efficient one-step solution
- Fast operating
- No preparation costs for sterilization, maintenance and service
- Low operating costs

Order number

Product

Description

ICT5011



**i-Cut**  
Power Morcellator for extracting tissue in LASH and myomektomy

i-Cut consists of:  
Power morcellator and obturator,  
1 piece, delivered sterile


ICT5020



**i-Cut Power Supply Type C – Euro**  
Reusable power supply unit for i-Cut ICT5011

1 piece, delivered non-sterile  
Additionally available:  
ICT5030 i-Cut Power Supply Type G  
ICT5040 i-Cut Power Supply Type I  
ICT5050 i-Cut Power Supply Type A  
ICT5060 i-Cut Power Supply Type C - Korea

MCS5111

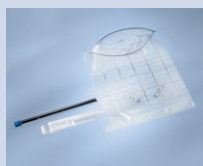

**More-Cell-Safe 0°**  
System for reducing the scattering of cancer cells during morcellation during laparoscopic surgery  
For 0° optics

The More-Cell-Safe system consists of:

- Cell-tight, intraperitoneally insufflatable More-Cell Bag with two accesses
- Visi-Shield to prevent camera contamination with 0° optics

10 systems per box, delivered sterile

MCS5151



**More-Cell-Safe 30°**  
System for reducing the scattering of cancer cells during morcellation during laparoscopic surgery  
For 30° optics

The More-Cell-Safe system consists of:

- Cell-tight, intraperitoneally insufflatable More-Cell Bag with two accesses
- Visi-Shield to prevent camera contamination with 30° optics

10 systems per box, delivered sterile

Proudly distributed by



☎ 1300 00 ENDO (3636)

info@endotherapeutics.com.au

www.endotherapeutics.com.au

Suite 301, 2 Banfield Road

Macquarie Park NSW 2113

**A.M.I. Headquarters:**

A.M.I. Agency for Medical

Innovations GmbH

Im Letten 1

6800 Feldkirch

Österreich

T +43 5522 90505 0

F +43 5522 90505 4006

e info@ami.at



www.ami.at

The products contained in this brochure meet the requirements of the Medical Device Directive 93/42/EEC and are accordingly marked with the CE mark.

CE 0297 ICT5011, MSC5111, MSC5151

EN 200052